


# Complete Sterile Wrap Protection

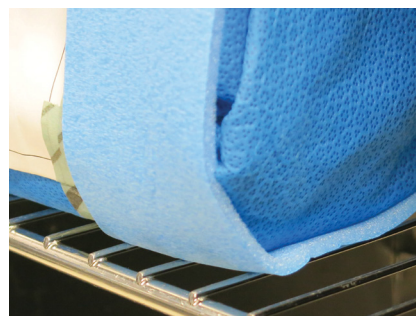

## Why use Tray Belts?

- Tray Belts, applied externally, protect the entire wrapped tray from damage caused by the environment during sterilization, storage and transport. Corner guards, placed inside the sterile wrap, only protect the sterile wrap from the tray's corners.
- May eliminate moisture problems, like wet packs, associated with padding placed within the sterile wrap.
- They prevent abrasion marks and damage caused by dragging wrapped trays, and by stacking wrapped trays on metal shelving units.
- They are compatible with steam sterilization.

Complete wrap protection at a lower cost.

The amount of Tray Belts material needed to protect a tray is relative to the tray size. Small and medium-sized trays cost less to protect. This isn't the case with corner guards, four corner guards are needed regardless of tray size. Large or small, there is no reduction in cost.

See page 2 for more detailed information.

SKU	DESCRIPTION	PER CASE
TBR3100NP	Tray Belts Roll, 3" x 100'	8 -100' Rolls and 1 roll of indicator tape
TBPC30	Tray Belts Pre-Cut Belts in a dispenser bag, 3" x 30"	200 Belts, (4 dispenser bags with 50 belts in each bag)
TBPC46	Tray Belts Pre-Cut Belts in a dispenser bag, 3" x 46"	200 Belts, (4 dispenser bags with 50 belts in each bag)
TBRD10	Tray Belts Acrylic Roll Dispenser – 7.5"w x 20"d x 11"h	1 Dispenser

Continued on the next page.

## CYGNUS MEDICAL


965 West Main Street, Branford, CT 06405 | Fax: 203.764.1858 | Email: sales@cygnusmedical.com  
www.cygnusmedical.com | 800.990.7489 **Patent Number 10751436**

## TBR3100NP Roll | Product Comparison – Tray Belts vs. Corner Guards

Does the Product...	Tray Belts	Corner Guards
Help to eliminate moisture collection problems (wet packs) that can compromise sterile wrapping during steam sterilization?	✓ Tray Belts are placed on the outside of the tray and do not collect moisture.	✗ <b>No!</b> Likely to cause wet packs.
Protect sterile wrapping from rips caused by the tray's corners?	✓	✓
Protect sterile wrapping from rips caused by stacking?	✓	✓ Only when the top tray has corner guards.
Protect sterile wrapping from abrasion marks caused by dragging trays off of sterilization racks?	✓	✗ No Protection
Protect sterile wrapping from sharp shelving edges?	✓	✗ No Protection
Protect sterile wrapping from rips caused by dragging trays out of case carts?	✓	✗ No Protection

## TBR3100NP Roll | Cost Comparison – Tray Belts vs. Corner Guards

SAMPLE TRAY SIZES				Tray Belts™ material needed to wrap 1 tray (with 2 belts.)*	Number of trays wrapped with 1 case of Tray Belts.**	COST PER TRAY		
						Tray Belts (2 Belts per tray)		Corner Guards (4 Guards per tray)
DIMENSIONS	W	L	H			@\$288.00 per case	@\$0.50 per guard	@\$0.60 per guard
<b>SMALL</b>	10"	10"	3"	36"	266	\$1.08	\$2.00	\$2.40
<b>MEDIUM</b>	10"	16"	3"	48"	200	\$1.44	\$2.00	\$2.40
<b>LARGE</b>	10"	20"	3"	56"	171	\$1.68	\$2.00	\$2.40
<b>EXTRA-LARGE</b>	10"	24"	3"	64"	150	\$1.92	\$2.00	\$2.40

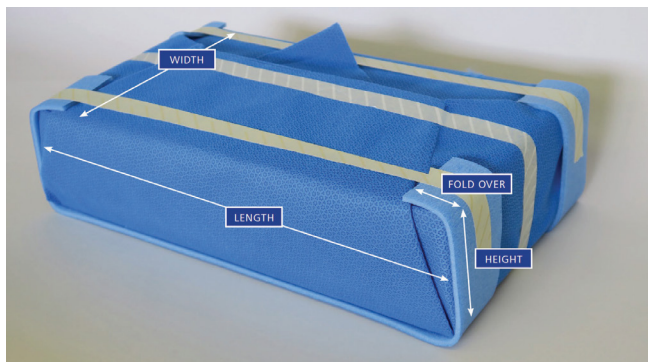

$$*(Length*2) + ((Height +1)*4)$$

\*\*A case of Tray Belts™ = 800 feet or 9600 linear inches of material.

### CYGNUS MEDICAL

965 West Main Street, Branford, CT 06405

www.cygnusmedical.com | 800.990.7489

**Patent Number 10751436**